

No Internet? No Problem.

Streamline your next wellness screening event with PTS Connect® ScreenPro, our portable data capture system designed specifically for wellness screening professionals. Unlike other systems, PTS Connect ScreenPro data hub does not require internet access and can deliver a customized wellness report for each participant using your choice of a company logo.



PTS Connect ScreenPro Can Help:

Save Money and Time

Run your events more smoothly and achieve faster throughput with PTS Connect ScreenPro. Our system does not require manual data entry of test results or internet access. Plug it in, press a button, and in a few moments, it's ready to go.

Safeguard Participants' Data

Give your customers and participants peace of mind with our secure, self-contained, multiple-encrypted, HIPAA-compliant data capture system.

Drive Better Wellness Decisions

Participants want to leave a screening event with results they can understand and can act upon. With PTS Connect ScreenPro, results are available for online review within moments, with the option of creating a printable report.

Powered by Oxygen


The PTS Connect ScreenPro data hub was developed specifically for the professional screening market using Oxygen software—a multiple-encrypted software and hardware solution for remote and non-WiFi enabled environments.

Give Your Customers What They Need Most

Companies and organizations want accurate and well-organized aggregate data reports so they can evaluate the efficacy of their wellness programs. These reports can be easily downloaded to a PC at the end of your screening event.

Get Started Today!

Don't waste another day or another dollar with outdated and cumbersome data management systems. Learn how you can get started with PTS Connect ScreenPro by calling +1 (317) 870-5610.

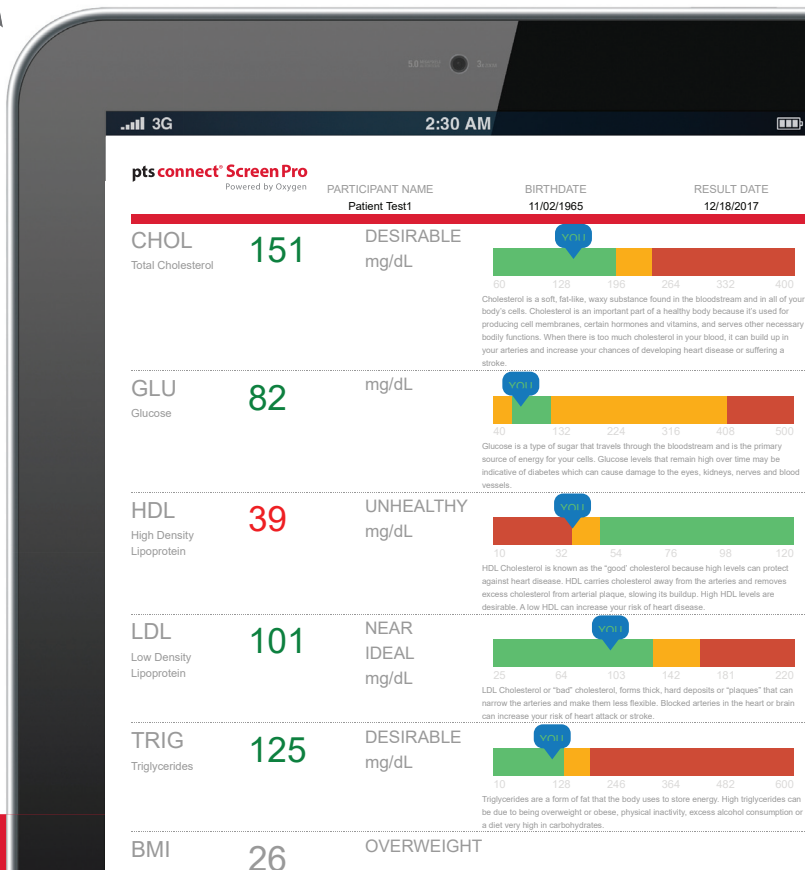




How It Works


- 1)** After powering on, the PTS Connect ScreenPro data hub creates a secure intranet that allows for wireless communication via the UBIQUITI® Access Point with tablets, PCs, CardioChek® Plus analyzers, and select printers.



- 3)** After participants are registered, they are automatically placed into a queue for testing. With our proprietary wireless integration, test results are automatically sent to the PTS Connect ScreenPro data hub where they are merged with the participants' records.

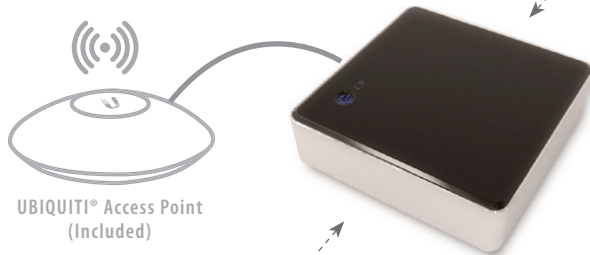

- 4)** After participants have completed the screening, they can review their results. If desired, the system can wirelessly print individual results onsite so that the participant may leave with a copy.

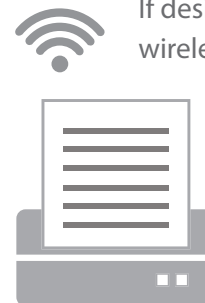

- 2)** Once the wireless network is ready, tablet computers are used to register participants. As data is entered, it is automatically transferred to the data hub for storage, ensuring no participant data remains on the tablet.



- 5)** After all participants have completed the screening event, reports can be easily generated. These reports can be downloaded to a PC, and then securely uploaded to the company servers or portal.



For more information about the PTS Connect ScreenPro, call +1(317) 870-5610 today.